


Imaging part

Image processing, display and recording

All-in-One



Specification	
Model	VLM-03
Display size	8 inches
Weight	1.0 Kg
Storage	32 GB
Display resolution	1024*768 pixels
Video output capacity	HDMI, 1920*1080 pixels
Battery capacity	≥6400 mAh
Battery working time	Up to 4 hours

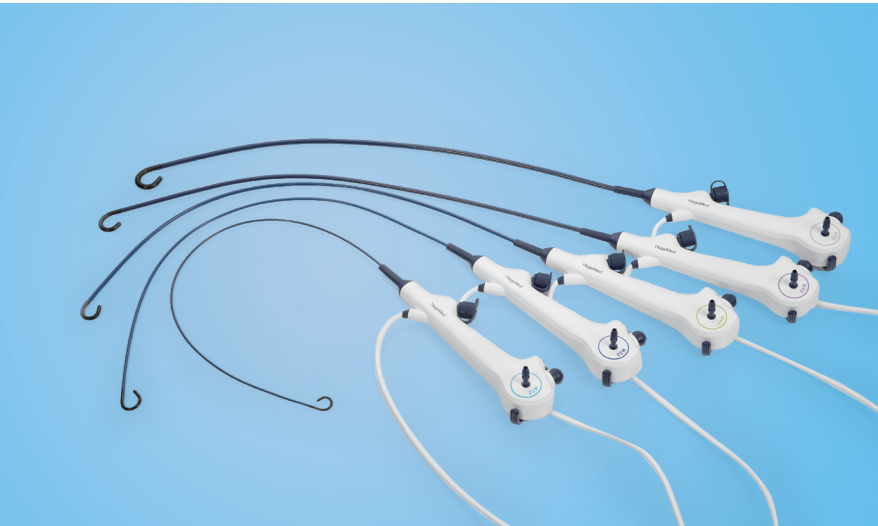

Brochure
**Single-use
Bronchoscope**



Ordering information

Product	Item code
Single-use bronchoscope	
BR-M22	2016B-PA005
BR-M32	2016B-PA004
BR-M40	2016B-PA003
BR-M50	2016B-PA001
BR-M58	2016B-PA006
Image processor	
VLM-03	2021B-PA003

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

HugeMed

Shenzhen HugeMed Medical Technical Development Co., Ltd.
415, 416-1, 516-1, Building 2, No.1, Mawu Road, Baoan Community, Yuanshan Street,
Longgang District, Shenzhen, Guangdong, 518115, China
+86-755-22275833 +86-755-22275833 info@hugemed.net

For more information
please visit our website:
www.hugemed.net



HugeMed

P/N: ENG-EU-CATA-fBRS-suBR-INTENT-202306V1.0

An optimal Bedside Bronchoscopy Solution



Portable

- Light weight and compact;
- Especially convenient for ICU daily examination.



Cost-effective

- Raises cash flow ratio;
- No maintenance cost;
- No disinfection cost.



No Cross-infection

- Single-use, one scope for one patient, no cross infection trouble.



Always Ready To Go

- No need to wait for reprocessing, always available for your needs.



More Sizes More Possibilities

HugeMed™ BR series single-use bronchoscope has 5 models that targets different clinical usages. By a feat of an ultra-compact CMOS sensor, the scopes are having slimmer outer sizes but larger working channels, which provides more possibility to finish different therapy operations.

Single-use Bronchoscope



HugeMed™ BR Series
BR-M58
5.8/3.4mm



HugeMed™ BR Series
BR-M50
5.0/2.8mm



HugeMed™ BR Series
BR-M40
4.0/2.0mm



HugeMed™ BR Series
BR-M32
3.1/1.2mm



HugeMed™ BR Series
BR-M22
2.2/0mm



Deflection	Working Length	Field of View
210° Up/Down	600mm	120°



Chip on tip
CMOS+LED



Ergonomic Design
Ease-of-use



Intensive Bending
210° up/down



Single-use
Sterilized

BR-M22
2.2/0mm

BR-M32
3.1/1.2mm

BR-M40
4.0/2.0mm

BR-M50
5.0/2.8mm

BR-M58
5.8/3.4mm

SLIM MODELS

- A petite person intubation or airway examination;
- DL tube position checking;
- Reach to smaller bronchi for lung peripheral checking or therapy.

ROUTINE MODELS

- The size of BR-M40 and BR-M50 can pass most of the adult ET tube, suitable for routine intubation and bedside bronchoscopy therapy. The thickness of their insertion tubes offers more maneuverability when doing the operation.

THERAPY MODELS

- The scopes whose working channels larger than 2.8mm normally are defined as the therapy models. HugeMed™ BR-M58 has a very fine working channel of 3.4mm, which provides 28.5% space more than that of 3.0mm, and 47.4% than that of 2.8mm.



BR-M50
5.0/2.8 mm

Compare with another flexible bronchoscopes whose outer diameter is nearby 5.0mm but the working channel is less than 2.6mm, the BR-M50 is a versatile bronchoscope with a regular outer diameter but a large channel which brings the advantages of:

- Able to reach smaller bronchus;
- Adapt to more ETT/DLT sizes;
- More power for suction.